

Cap-Chun

Table of Contents

Alphabetized Listing

A		
Arrow Adaptor		47
B		
Blank Adaptor.22		47
C		
Carrying Case - Pistol		47
Carrying Case - Rifle		48
Charges		50
CO2 Gas Packs		51
G		
Gas Packs		51
Guns and Accessories		47-48
M		
Marking Syringe Complete		50
N		
Needles & Syringes		48-50
Noseplug Barb		48
Noseplug Heavy		49
Noseplug Large Collar		49
Noseplug Plain		49
Noseplug Small Collar		49
Noseplugs		48-49
P		
Pistol		47
Plunger Lube		48
Positioner		48
Powerloads		51
Powerloads - Low		51
Powerloads - Medium		51
Powerloads - Very Low		51
Powerloads & Charges		50-51
Practice Syringe		50
Projectors		47
Projectors - Accessories		47-48
R		
Rifle - Long Range		47
Rifle - X-Long Range		47
Rubber Plunger		48
S		
Sideport		49
Sideport Noseplug - Drop Out		49
Sideport Small Collar		49
Stock Extension		48
Syringe Barrel Only		50
Syringes		50
T		
Tail Piece		50



Guns and Accessories

Projectors

Pistol



Short range projector: a pistol type projector which is powered by one CO2 gas pack and can fire approx. 20 syringes for a distance of 40 ft. each. Mid range projector: a pistol type projector that is breech loaded and is powered by 2 - 12 gram CO2 gas packs.

KVS	Size	Mfg.
039-305	Short Range	Cap-Chur
039-307	Mid Range	Cap-Chur

Rifle - Long Range



Long range projector: a rifle type projector which is powered by two CO2 gas packs. The long range projector can be expected to fire 12 x 1cc syringes a distance of 25-30 yards with improved accuracy over the short range pistol.

KVS	Size	Mfg.
039-304	Long Range	Cap-Chur

Rifle - X-Long Range



This rifle type projector is breech loaded and is powered by a .22 cal. Blank cartridge. It comes with Blank adaptor for the .22 cal. Power loads. Multiple ranges are achieved through the use of 4 different Power Loads. This model has a .50 cal., rifled barrel. The Barrel length is 25 inches, the overall length is 41 inches and the weight is 7lbs. It can be fitted with a scope. Used for ranges of 5-75 yards.

KVS	Size	Mfg.
039-302	X-Long Range	Cap-Chur

Projectors - Accessories

Arrow Adaptor



This adaptor enables you to attach a Cap-Chur syringe to an arrow, for use with cross or longbows.

KVS	Mfg.
039-465	Cap-Chur

Blank Adaptor.22



The Blank Adaptor is used with the extra-long range projector. It holds the power-load inside the breech.

KVS	Mfg.
039-449	Cap-Chur

Carrying Case - Pistol



Vinyl-covered, padded, carrying case. Zippered closure. Convenient pockets for storing your syringes and equipment.

KVS	Size	Mfg.
039-463	Pistol	Cap-Chur



Guns and Accessories

Projectors - Accessories

Carrying Case - Rifle



Vinyl-covered, padded, carrying case. Zippered closure. Convenient pockets for storing your syringes and equipment.

KVS

039-461

Size

Rifle

Mfg.

Cap-Chur

Plunger Lube



Premium silicone grease used to allow the rubber plunger to move forward freely and seal tightly along the barrel wall.

KVS

039-321

Size

1.5oz

Mfg.

Cap-Chur

Positioner



POSITIONER: is a clear fibreglass rod used to position the rubber plunger inside the syringe barrel.

KVS

039-445

Mfg.

Cap-Chur

Rubber Plunger



RUBBER PLUNGER: Comparable to the plunger on any syringe. Seals against the inside barrel wall and forces medication out through the needle.

KVS

039-443

Mfg.

Cap-Chur

Stock Extension



Shoulder stock for short range pistol.

KVS

039-309

Mfg.

Cap-Chur

Needles & Syringes

Noseplugs

Noseplug Barb



NC NEEDLES: used primarily on thin skinned animals that are being "knocked out". Identified by the barb on the needle shank. Needle acts like a fish hook and will not come out in spite of animals activity.

KVS

039-383

Size

3/4"

Part

NC1

Mfg.

Cap-Chur

039-385

1 1/8"

NC2

Cap-Chur

039-387

1 1/2"

NC3

Cap-Chur

039-388

1 3/4"

NC4

Cap-Chur

039-390

2"

NC5

Cap-Chur

039-391

2.5"

NC6

Cap-Chur



Needles & Syringes

Noseplugs

Noseplug Heavy



Noseplugs (needles) for thick skinned animals. (500lbs & over). Diameter of the needle measures 0.125".

KVS	Size	Part	Mfg.
039-401	1 1/4"	NCL1	Cap-Chur
039-403	1 1/2"	NCL2	Cap-Chur
039-405	1 3/4"	NCL3	Cap-Chur

Noseplug Large Collar



NS NEEDLE: Used for subcutaneous injections, identified by the short needle shank.

KVS	Size	Part	Mfg.
039-371	5/8"	NC1	Cap-Chur
039-373	3/4"	NS2	Cap-Chur

Noseplug Plain



PLAIN NEEDLES: Smooth shanked (no barb or collar) needle.

KVS	Size	Mfg.
039-389	3/4"	Cap-Chur
039-393	1 1/8"	Cap-Chur

Noseplug Small Collar



NM NEEDLES: Used for IM injections, identified by the longer needle shank. The collar near the end of the NS & NM needles helps ensure the needle does not "pop out" of the animal before the injection is complete.

KVS	Size	Part	Mfg.
039-375	1 1/8"	NM1	Cap-Chur
039-377	1 1/4"	NM2	Cap-Chur
039-379	1 1/2"	NM3	Cap-Chur
039-381	1 3/4"	NM4	Cap-Chur

Sideport

Sideport Noseplug - Drop Out



Sideport needle with drop out design.

KVS	Size	Mfg.
039-430	3/4"	Cap-Chur

Sideport Small Collar



NEW: The sideport delivery needle is designed primarily to eliminate the problem of "bounce outs" due to back pressure. Medication is dispensed over an area of about four times that of a regular needle.

KVS	Size	Part	Mfg.
039-425	1 1/8"	NM1	Cap-Chur
039-427	1 1/4"	NM2	Cap-Chur
039-429	1 1/2"	NM3	Cap-Chur



Needles & Syringes

Syringes

Marking Syringe Complete



For marking animals with colored ink using the Cap-Chur projectors as the delivery system. Marking syringes have a flat face noseplug, which does not pierce the skin of an animal. Upon impact with the animal, ink is squirted out of the exhaust ports in the syringe barrel walls, colouring the animals hide.

KVS	Size	Mfg.
039-475	15ml	Cap-Chur

Practice Syringe



Practice syringe for the Cap-chur gun.

KVS	Size	Mfg.
039-447	1-3ml	Cap-Chur

Syringe Barrel Only



The Cap-Chur syringe barrel is the holding tank or carrier of the injection while the syringe is in flight. The diameter of the syringe barrel must remain constant in order to fit through the projector. Therefore the length of the syringe barrel is proportional to how many ml's the barrel holds.

KVS	Size	Mfg.
039-347	1ml	Cap-Chur
039-349	2ml	Cap-Chur
039-351	3ml	Cap-Chur
039-353	4ml	Cap-Chur
039-355	5ml	Cap-Chur
039-357	7ml	Cap-Chur
039-359	10ml	Cap-Chur
039-361	15ml	Cap-Chur
039-363	20ml	Cap-Chur

Tail Piece



TAIL PIECE: is the "feather" or guide for the Cap-Chur syringe. Helps keep the syringe "true" in flight.

KVS	Mfg.
039-441	Cap-Chur

Powerloads & Charges

Charges



Charges are very small explosive devices which "fire" on impact. When the syringe hits the target, the rubber plunger is pushed forward by the expanding gas caused by the Cap Chur charge explosion. The drug is then forced out through the needle and into the target.

KVS	Size	Pack	Mfg.	TDG Class
039-451	1-3ml	50	Cap-Chur	LQ
039-453	4-10ml	50	Cap-Chur	LQ
039-455	15ml	50	Cap-Chur	LQ

CAPCHUR



Powerloads & Charges

Gas Packs

CO2 Gas Packs



Used in the long range + short range projectors to propel the Cap Chur syringes.

KVS
156-001

Pack
5

Mfg.
KVS

Powerloads

Powerloads - Low



POWER LOADS: .22 calibre blank loads used to propel Cap Chur syringes out of the Cap Chur Extra Long Range Projector. LOW: used at ranges from 15-40 yds.

KVS
039-313

Color
Green

Size
Low

Pack
50

Mfg.
Cap-Chur

TDG Class
LQ

Powerloads - Medium



POWER LOADS: .22 calibre blank loads used to propel Cap Chur syringes out of the Cap Chur Extra Long Range Projector. MEDIUM: used at ranges from 30-60 yds.

KVS
039-315

Color
Yellow


Size
Medium

Pack
50

Mfg.
Cap-Chur

TDG Class
LQ

Powerloads - Very Low

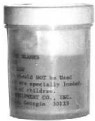

POWER LOADS: .22 calibre blank loads used to propel Cap Chur syringes out of the Cap Chur Extra Long Range Projector. VERY LOW: used for ranges from 5-25 yds.

KVS
039-311

Color
Brown

Size
Very Low

Pack
50

Mfg.
Cap-Chur

TDG Class
LQ